



INNOVA
PRO

**IONIC
PRO C10**

INNOVA PRO
HEALTHCARE SCIENCE

IONIC PRO C10

INNOVA PRO SYSTEM TECHNOLOGY

IONIC PRO C10 is the most advanced diathermic device of the moment. Our IONIC system upgrades cutting edge technology with a PFC power control to perform the best treatment results.

Our device develops 3 different unique frequencies for multiple applications. We are a market reference due to the variety of applicators available and treatment possibilities.

IONIC PRO C10 BY INNOVAPRO
UNCOMPARABLE DIATHERMY DEVICE
DISCOVER IT...



WHAT IS A DIATHERMY?


Diathermic treatments are for beauty, wellness, and health.

The etymological definition is = DIA (through or via) + TERMIA (heat)

The principle consists of the physiological activation of the deep tissues by 3 phenomena: the micro circulation, vessel dilatation and increase of local temperature.

Through heating we obtain a stimulation on the cell activity that will promote the generation of collagen and elastin, an improvement in the blood flow, an increase oxygenation that will improve and rejuvenate the skin. By the same way we can achieve a decrease of the body fat and a lifting effect on the skin.

Diathermy is the ideal treatment to rejuvenate, firm and reshape the skin.

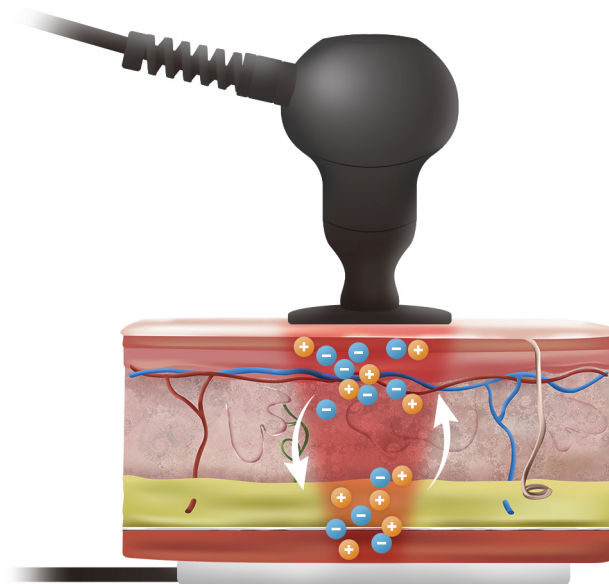

The general effects of diathermy are as follows:

- ✓ Activation and restructuring on the production of collagen and elastin.
- ✓ Improvement on skin regeneration and in terms of hydration, elasticity, and smoothness.
- ✓ Rejuvenation and revitalization of the skin.
- ✓ Body fat removal. Cellulite removal.
- ✓ Improvement of vascularization.
- ✓ Healing of muscle contracture and articular mobility.
- ✓ Reduction of edema and local pain.



OUR DEVICE IONIC PRO C10 3 FREQUENCIES

Frequency of 469 KHz

IONIC PRO performs in the superficial layer with the CAPACITIVE probe and in the Deep layer with the RESISTIVE probe. Both work with and output power of 300W.



What effects do we get?

CAPACITIVE:

- Firms and improves the skin.
- Reduces wrinkles.
- Activates circulation.
- Oxygenates tissues.
- Eliminates cellulite an "Orange peel skin."
- Improves stretch marks.

RESISTIVE:

- Eliminates body fat.
- De-contractures muscles.
- Muscle stimulation.
- Firming of deep tissues.

Our device is the only one that can perform with 2 simultaneous capacitive applicators DUAL CONNECT C10



Frequency of 300 KHz.

IONIC PRO has the ability to perform with a 300 KHz frequency for Deep pain therapeutic treatments.

IONIC PRO also has a special applicator "The Horseshoe" ideal for our therapeutic, anti-cellulite and modeling treatments.



WHEN FREQUENCY IS THE KEY...



BIPOLAR 500 kHz

IONIC PRO C10 is equipped with a bipolar applicator with a 500 KHz frequency output combined with chromo therapy. It is the strongest and deepest applicator in the market. It is the ideal complement to achieve an intense flash effect on our rejuvenation treatments. Your patients will see the difference and so will you.

WHEN LOOKING GOOD IS IMPORTANT...



MAXIMUM OUTPUT OF ENERGY 300W (max)

- STABLE ENERGY OUTPUT.
- PEAK FREE OUTPUT.
- DEEP AND INTENSE RF TECHNOLOGY.
- CAPACITIVE, RESISTIVE AND MULTIPOLAR.
- FOR FACE AND BODY.
- IDEAL FREQUENCY OF 469 KHz.
- MULTIPLE ACCESORIES.

CAP



RES



Bipolar
with cable



Dual
connector



Return plate



Return bar



STIMULATES
RESHAPES
REESTRUCTURES



REF: SNE-7100

TECHNICAL SPECIFICATIONS IONIC PRO C10

Diathermy with Capacitive and Resistive transfer.

Voltage:	100-250 VAC-50/60Hz
Source:	350 VA +- 10%, 320W, 1.265°
Output frequency:	CAP Deep: 300 KHz and CAP 469KHz and RES 500 KHz (0.3 MHz, 0.4MHz y 0.5 MHz) +-10%
Channels:	3 (1 -CAP, 2 -RES, 3 -BIPOLAR) + Return cathode
Energy Output	RES 300W (50Ω de impedance) mode A RES 210W (50Ω de impedance) mode B RES 255W (50Ω de impedance) mode C CAP 320 W (200Ω de impedance) normal mode (650 VA) CAP 380 W (200Ω de impedance) Deep mode (655 VA) BIPOLAR 130W (10Ω de impedance)
Power Factor Controller (P.F.C)	0.99 in oscilloscope
Accessories	Monopolar (Resistive RET), Monopolar (Capacitive CAP), Bipolar, Capacitive Dual, Horseshoe, return cathode, return bar.
Applicator size	Capacitive: 20-30-40-60 -80mm (+capacitive doble 60mm) Resistive: 30-60-80mm (+protector 60mm)
Screen	Tactile 7" (800xRGBx480)
Protection level	Class I ME, type BF
Weight	13, 18 Kg
Size	400x404x159 mm
Certification	CE2460 marking according to MDD93/42/EEC class IIa, RoHS, WEEE, ISO13485





TOGETHER
WE CREATE BEAUTY

CENTRAL HEADQUARTERS

Street 16
Warehouse 39-41-43
Industrial Park Base 2000
Lorquí - Murcia

968 603 256

www.innovapro.es

info@innovapro.es

 [innovapro_oficial](https://www.instagram.com/innovapro_oficial)  [InnovaPro](https://www.facebook.com/InnovaPro)